



# MINI Water flosser




Product Manual



## MINI Water flosser

This compact, travel-ready mini water flosser is designed to give you a deep clean that's gentle on gums but tough on plaque. With a lightweight and portable design, USB charging, and just the right tank size for efficient flossing, it's the perfect companion for life on the go.

### SAFETY CAUTION

 NOTE: PLEASE READ THIS INSTRUCTION MANUAL CAREFULLY BEFORE USE.

IT CONTAINS IMPORTANT OPERATING INFORMATION AND SAFETY PRECAUTIONS. IF YOU HAVE ANY EXISTING ORAL HEALTH CONDITIONS OR CONCERNS, SEEK ADVICE FROM YOUR DENTAL OR MEDICAL PROFESSIONAL PRIOR TO USE.

### WHAT'S INCLUDED


- Mini Water Flosser (Flosser device, tank, and tube)
- USB charging cable (Only use the provided USB cable for charging product.)
- 4 nozzle attachments
- User manual

### DISCLAIMER

This product is not intended to diagnose, treat, cure, or prevent any medical or oral health condition.

If you have concerns about your oral health, please consult a qualified dental or medical professional before use.

### **MINI WATER FLOSSER**

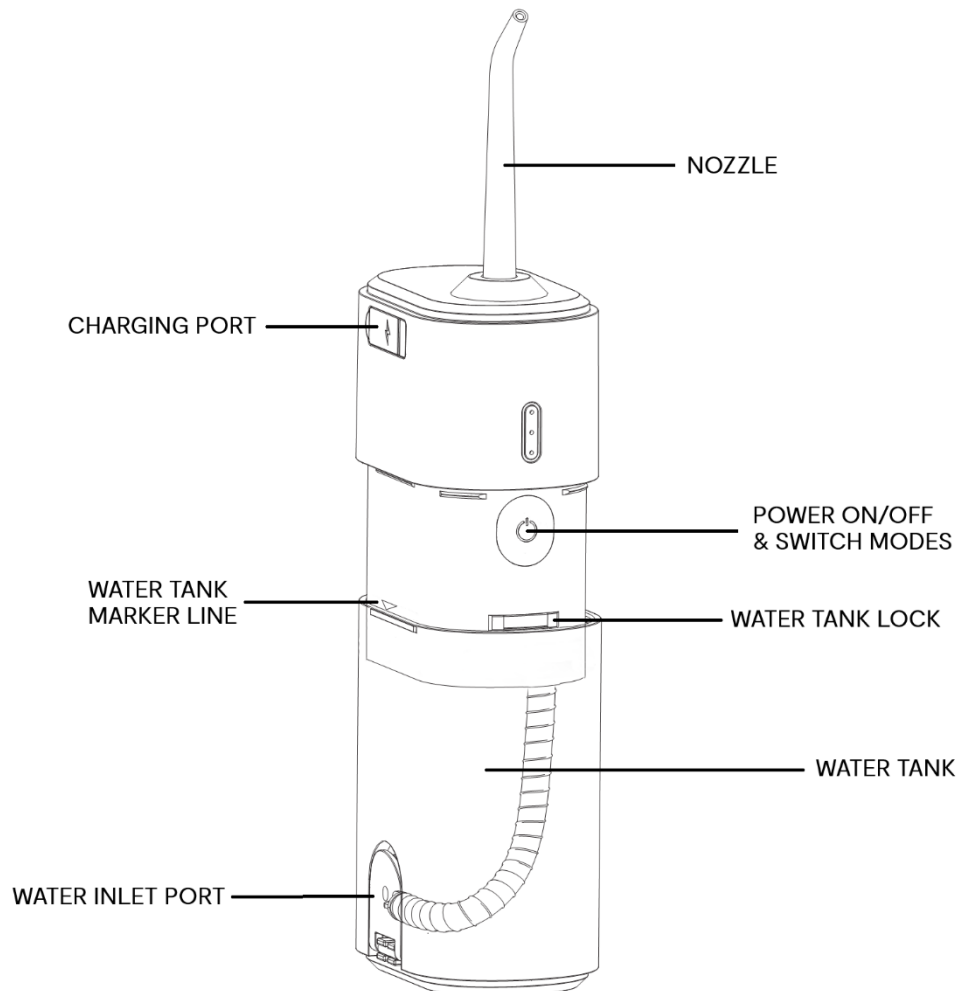

### **IMPORTANT INFORMATION BEFORE YOU BEGIN**

- Keep out of reach of children.
- Use only cold or warm water in the tank.
- Always fill the water tank before turning the device on.
- Do not operate the device with an empty water

tank.

- Stop using the device immediately if it becomes damaged.
- Keep the charging cable away from heated surfaces.
- Do not add salt, iodine, bleaching agents, or insoluble essential oils to the tank, as this may affect performance or shorten the product life.
- This device contains a rechargeable lithium battery.
- If you notice unusual or persistent gum bleeding while flossing, stop use and consult your dentist.
- If you are currently receiving treatment for any oral health condition, seek advice from your dental professional before using this product.

#### **NOZZLE TIPS**

THE CLASSIC TIP



For everyday use. It cleans deep between teeth and below the

PIK POCKET TIP



A soft rubber tip designed for gentle, low-pressure rinses below the

THE PLAQUE SEEKER TIP



With 3 thin tufts of bristles, this tip helps to remove plaque around hard-to-reach areas.

TONGUE CLEANER TIP



A spoon-like shape that helps to clean

#### **HOW TO CHARGE**

1. Insert the USB cable into a 5V USB charging port and

- connect it to the device.
2. Three blue lights will begin flashing, indicating that the device is charging.
  3. Once fully charged, the three lights will stay on without flashing.

## **PRESSURE MODE**

### A) WEAK

A gentle, low-pressure mode ideal for beginners or sensitive gums.

We recommend starting on this mode.

### B) STRONG

A high-pressure mode for a deeper clean. Suitable for everyday use once you're comfortable.

### C) PULSE

A rhythmic pulsing mode designed to massage and stimulate the gums.

## **HOW TO USE**

**IMPORTANT:** Always ensure the water tank contains water before turning the device on. Do not operate the device with an empty tank, as this may cause pressure build-up and make the tank more difficult to remove. Before use, fully extend the water tank past the water tank marker line and confirm it is securely locked in place by the water tank lock.

### **1. Open the air-release flap**

Before opening or closing the water tank, open the water inlet valve at the bottom of the tank. This allows air to enter the tank and prevents water inlet valve vacuum pressure.

### **2. Fully extend the water tank**

Hold the top of the device firmly and slide the water tank downward until it is fully extended. Ensure the water tank is positioned past the water tank marker line and secured in place by the water tank lock before use.

### **3. Fill the tank with water**

Fill the tank with cold or warm water only. Do not overfill.

**4. Close the air-release flap**

Once filled, close the water inlet valve securely before use.

**5. Attach a nozzle tip**

Insert the desired nozzle firmly into the nozzle port.

**6. Position the device**

Place the nozzle tip in your mouth and aim along the gumline.

**7. Turn the device on**

Press the power button once to start. Do not turn the device on unless the tank contains water.

**8. Select pressure mode**

Press and hold the power button for 2 seconds to switch between pressure modes.

**9. Floss correctly**

Slowly guide the nozzle along the gumline and between teeth. Keep your mouth slightly open during use to avoid splashing.

**10. Maintain correct position**

Keep the device upright at all times to maintain proper water flow.

**11. Refill if needed**

Refill the water tank as needed to complete your flossing session. Always turn the device off before refilling.

**12. Turn off after use**

Press the power button to turn the device off after use. Do not operate the device with an empty water tank.

**CLEANING AND MAINTENANCE**

A) Nozzle

Rinse with clean water and dry with a soft cloth or paper towel.

B) Device

Ensure the power is turned off. Disconnect the charger and seal the charging port. Open the water tank, empty its contents, and close it again.

Turn the device on briefly to release any remaining water, then turn it off.

Wipe the device dry with a soft cloth.

C) Water tank

We recommend cleaning the water tank once a week. After each use, empty any remaining water from the tank, open the water inlet valve, and allow the water flosser to air dry completely before storage. This helps reduce moisture build-up inside the tank. For cleaning, rinse the tank with clean, warm water only.

#### **CAUTION**

- This product is not waterproof. Do not submerge the device in water.
- Keep the device away from children, direct sunlight, high temperatures, and humidity.
- Consult your healthcare professional before use if you have any underlying health conditions.
- Do not charge the device while water is in the tank. Always unplug the charger before filling the tank, and do not overfill.
- Do not operate the device without water in the tank, as this may cause damage.
- Do not operate the device without a nozzle tip attached.
- Do not aim or spray water directly under the tongue, or into the eyes, ears, nose, or other parts of the body.
- Do not use this product if you have open wounds or sores in your mouth.
- Remove all jewellery inserts from your mouth before use.
- Use this product only for oral cleaning purposes.
- Do not operate the device without a nozzle tip

attached.

#### **PRODUCT SPECIFICATIONS**

Product name	OXX Bodycare Mini Water Flosser
Model no.	NTS.802
Output voltage	3.7V
Input voltage	5V



#### **LIABILITY**

The sponsor and/or manufacturer does not accept responsibility for the misuse of this product. It is the customer's responsibility to read and understand all instructions and safety information in this manual and to follow all guidelines carefully. The sponsor and/or manufacturer will not be liable for any personal injury, property damage, or death resulting from improper use of the device.