


SIZE:44.1X21cm

anko

Heat Pad Small

Instruction Manual

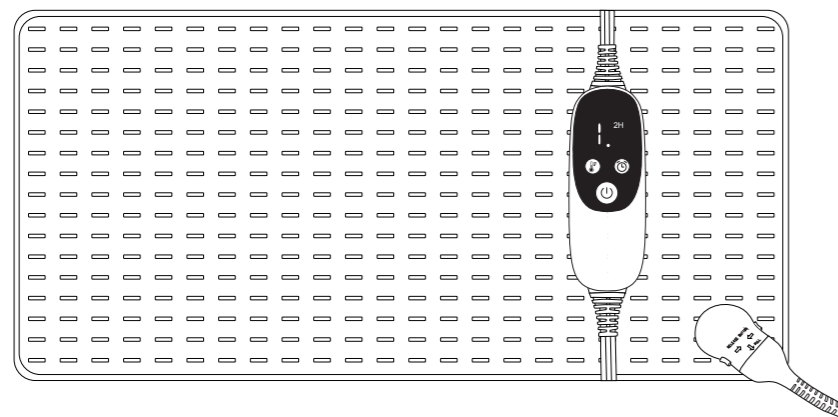

Model H-5525-5

To ensure proper operation of this device, please read this instruction manual carefully before using the device and keep it for further reference.

(1.)

READ AND SAVE IMPORTANT SAFETY INSTRUCTIONS

DANGER- TO REDUCE RISKS OF BURNS, ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

1. DO NOT USE WHILE SLEEPING. DO NOT USE WHILE BATHING OR IN A SHOWER.
2. DO NOT USE ON AN INFANT OR ON AN ANIMAL.
3. CHILDREN UNDER THE AGE OF THREE ARE NOT TO USE THIS APPLIANCE DUE TO THEIR INABILITY TO REACT TO OVERHEATING.
4. THE APPLIANCE CAN BE USE BY CHILDREN OLDER THAN 3 YEARS AND YOUNGER THAN 8 YEARS UNDER SUPERVISION AND WITH THE CONTROL ALWAYS SET TO MINIMUM TEMPERATURE VALUE.
5. THE APPLIANCE IS NOT TO BE USE BY YOUNG CHILDREN OVER THE AGE OF THREE UNLESS THE CONTROLS HAVE BEEN PRE-SET BY A PARENT OR GUARDIAN, AND UNLESS THE CHILD HAS BEEN ADEQUATELY INSTRUCTED ON HOW TO OPERATE THE CONTROLS SAFELY.
6. THIS APPLIANCE CAN BE USE BY CHILDREN AGED FROM 8 YEARS AND ABOVE AND PERONS WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES OR LACK OF EXPERIENCE AND KNOWLEDGE IF THEY HAVE BEEN GIVEN SUPERVISION OR INSTRUCTION CONCERNING USE OF THE APPLIANCE IN A SAFE WAY AND UNDERSTAND THE HAZARDS INVOLVED.
7. CHILDREN SHALL NOT PLAY WITH THE APPLIANCE. CLEANING AND USER MAINTENANCE SHALL NOT BE MADE BY CHILDREN WITHOUT SUPERVISION.
8. THIS PAD IS NOT TO BE USED ON OR BY AN INVALID, A SLEEPING OR UNCONSCIOUS PERSON, A PERSON WITH DIABETES, OR A PERSON WITH POOR BLOOD CIRCULATION.
9. THIS PAD MUST NOT BE USED BY PERSONS INSENSITIVE TO HEAT AND OTHER VERY VULNERABLE PERSONS WHO ARE UNABLE TO REACT TO OVERHEATING.
10. IF THE APPLIANCE IS SLEPT ON WITH THE CONTROLS SET TO A HIGHER TEMPERATURE THE USER MAY SUFFER SKIN BURNS OR HEAT STROKE.
11. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING. CHECK SKIN UNDER PAD FREQUENTLY. DO NOT USE IN AN AREA WHERE OXYGEN IS BEEN ADMINISTERED.
12. DO NOT INSERT PINS INTO THE ELECTRIC PAD, THEY MAY DAMAGE THE ELECTRIC WIRING.




(2.)

13. USE HIGH SETTING TO PRE HEAT BEFORE USE ONLY. DO NOT USE WITH THE CONTROL SET TO THE HIGH SETTING, IT IS HIGHLY RECOMMENDED THAT THE PAD TO BE SET TO LOW HEAT FOR CONTINUOUS USE.
14. THE CORD AND CONTROL MUST BE AWAY FROM HEAT SOURCES SUCH AS HEATING AND LAMPS.
15. DO NOT BUNCH, SIT ON OR AGAINST OR CRUSH PAD, AVOID SHARP FOLDS.
16. EXAMINE THE APPLIANCE FREQUENTLY FOR SIGNS OF WEAR OR DAMAGE. IF THERE ARE SUCH SIGNS, IF THE APPLIANCE HAS BEEN MISUSED OR DOES NOT WORK, DO NOT USE IT.
17. THIS APPLIANCE IS NOT INTENDED FOR MEDICAL USE IN HOSPITALS.
18. PLACE PAD ON TOP OF AND NOT UNDER THE PART OF THE BODY NEEDING HEAT, DO NOT LAY ON THE PAD. DO NOT USE THE PAD FOLDED.
19. USE THIS PAD ONLY AT RATED VOLTAGE, AS MARKED ON THE PRODUCT, DO NOT USE WITH ELECTRIC POWER GENERATORS, POWER CONVERTERS OR INVERTERS. UNPLUG WHEN NOT IN USE.
20. DO NOT TAMPER WITH THIS PAD IN ANY WAY.
21. THERE ARE NO USER-SERVICEABLE PARTS IN THE PAD.
22. TO AVOID SKIN BURNS, DO NOT USE THIS PAD WITH LINIMENTS, SALVES OR OINTMENT PREPARATIONS, ESPECIALLY ONES THAT CONTAIN HEAT-PRODUCING INGREDIENTS. SKIN BURNS COULD RESULT.
23. THE CABLE AND CONTROL UNIT OF THE APPLIANCE CAN GIVE RISE TO RISKS OF ENTANGLEMENT, STRANGULATION, TRIPPING OR TREADING IF NOT CORRECTLY ARRANGED. THE USER SHALL MAKE SURE THAT EXCESS TIES AND ELECTRIC CORDS SHALL BE ARRANGED IN A SAFE WAY.
24. IF THE SUPPLY CORD IS DAMAGED, IT MUST BE REPLACED BY THE MANUFACTURER, ITS SERVICE AGENT OR SIMILARLY QUALIFIED PERSONS IN ORDER TO AVOID A HAZARD.
25. DO NOT USE THE PAD IF IT IS WET OR MOIST.
26. DO NOT OPERATE UNDER BLANKET OR PILLOW, EXCESSIVE HEATING CAN OCCUR AND CAUSE FIRE, ELECTRIC SHOCK, OR INJURY TO PERSONS.
27. NEVER OPERATE THIS PAD IF IT HAS A DAMAGED CORD OR PLUG, IF IT IS NOT WORKING PREPERLY. IF IT HAS BEEN DROPPED OR DAMAGED OR DROPPED INTO WATER.
28. NEVER LEAVE THE APPLIANCE UNATTENDED, ESPECIALLY WHERE CHILDREN ARE PRESENT.
29. LOOP CORD LOOSELY WHEN STORING. TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS.
30. WHEN STORING THE APPLIANCE, ALLOW IT TO COOL DOWN BEFORE FOLDING, DO NOT CREASE IT BY PLACING ITEMS ON TOP OF IT, STORE IT IN THE ORIGINAL PACKAGING, IN A DARK, COOL, DRY PLACE.
31. THE APPLIANCE IS ONLY TO BE USED WITH THE TYPES THAT ARE MARKED ON THE APPLIANCE.


(3.)

HOW TO USE YOUR HEATING PAD:

Make sure Control Connector is securely attached to the Pad Connector. Push both parts together for a tight fit. NOTE: control connector must be connected to heating pad before plugging unit into a power outlet socket.

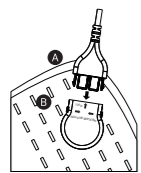

A. Control Connector

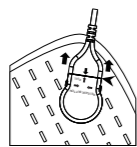


B. Pad Connector



C. Plug into the wall outlet

- After control connector is connected to the heat pad, plug the supply cord plug into a mains power outlet socket.
- To turn the pad on, press the **Power Button** (⏻). The screen will show "0" **NO HEAT**.
- Press the **HEAT SETTING** (🌡️) button to select the desired heat setting, 1 is the **low** heat, 2 is **medium** heat, 3 is **high** heat.
- Timer Setting:**
Press the **Timer Button** (🕒) to select the desired operating hours. Timer by default is set at **1H** when Heat setting on, press once for timer setting **1.5H**, press again for timer setting **2H**, then press to cycle back around to timer **1H**, Repeatedly press to increase the timer.
The numbers **1-1.5-2** are shown in hours and H represents hours.
The product will automatically shut off when the timer is up.
The heat setting and timer will loop when pressing the relevant function button.
- To turn the pad **OFF**, press the **Power Button** (⏻) one more time, the heating pad and controller display will turn off.
"3" is the highest setting and not recommended for prolonged use, we recommend using "low heat" setting for continuous use of the heat pad.
- DO NOT SIT ON OR AGAINST OR CRUSH THE PAD. AVOID FOLDING DURING USE.**
- DO NOT USE WHILE SLEEPING.**
- SEE ADDITIONAL WARNINGS.**

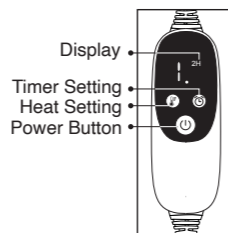

TECHNICAL SPECIFICATIONS:

Size: 61cm x 31cm
220-240V~ 50Hz 80W

CLEANING AND MAINTENANCE:

- Always turn OFF and unplug the heating pad for care and maintenance.
- Spot clean only.
- Do not use corrosive cleaning agents such as benzene, thinners, or solvents.
- Do not wipe the heating pad with a wet cloth - use a dry cloth only.

(4.)

- Do not immerse the product in water or any other liquid. If water is accidentally spilled on the product, it may cause electric shock.
- If the product will not be used for an extended period, cover it with a dust cover, remove the plug from the power outlet, and store it out of reach of children.
- Do not dispose of this product with household waste. Follow local environmental disposal regulations.

TROUBLE SHOOTING

Problem	Possible Reasons	Solutions
Pad doesn't get hot enough	The default heat setting is 0, the user did not select desired heat setting, heat level 3 is highest.	Use a desired heat setting, the heat setting can be adjusted from level 1~3.
Pad has a hot spot	The pad was folded or bunched up for a long time during use.	Make sure the pad is not folded or bunched up in any area.
The pad is turning yellow in spots.	Pad has a hot spot due to prolonged use of the HIGH setting with the heat pad folded or bunched.	Make sure the pad is not folded or bunched up in any area.
Pad has no heat	Check whether the controller and the pad are connected	1. Verify no physical damage to the cord, controller or pad. 2. Verify AC outlet is functional. 3. Plug the unit back into the AC outlet. 4. If the controller and the pad are well connected, it means the unit has been damaged, please do not use it again.
The display displays an "F" and blinks	1. Check whether the controller and the pad are connected. 2. If you plug in outlet before connecting the controller to the pad	1. Ensure the controller is connected to the pad before plugging into outlet. 2. If "F" appears, unplug the unit and try again following the correct connection sequence 3. If "F" still appears, that means the product has been damaged and cannot be used again.
The display appears "H" and blinks	High temperature protection	Disconnect the heat pad from mains power and allow it to cool down for 30 minutes, and then connect it back into the mains power outlet socket. Press the power button to restart.
The LED displays an " _ "	The time is up	Press the power button again to start

(5.)

12 Month Warranty

Thank you for your purchase.

Your new product is warranted to be free from defects in materials and workmanship for the period stated above, from the date of purchase, provided that the product is used in accordance with accompanying recommendations or instructions where provided. This warranty is in addition to your rights under the Australian Consumer Law.

For New Zealand customers, this warranty is in addition to statutory rights observed under New Zealand legislation.

We will provide you with your choice of a refund, repair (where possible) or exchange (availability dependent) for this product if it becomes defective within the warranty period. The business will bear the reasonable expense of claiming the warranty.

This warranty will no longer apply where the defect is a result of alteration, accident, misuse, abuse or neglect.

Please retain your receipt as proof of purchase and contact our Customer Service Centre as listed below, for the entity you purchased this product from, for any difficulties with your product. Warranty claims and claims for expense incurred in returning this product can be addressed to the respective Customer Service Centre.

Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and compensation for any other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure.

Contact for Kmart purchases	Contact for Target purchases
Kmart Australia Ltd. C/- Customer Service Centre L2/ 1 MIDDLE ROAD, CHADSTONE VIC 3148 AUSTRALIA Customer Service: 1800 124 125 (Australia) or 0800 945 995 (New Zealand) or via Customer Help at kmart.com.au	Target Australia Pty Ltd C/- Customer Service Centre L2/ 1 MIDDLE ROAD, CHADSTONE VIC 3148 AUSTRALIA Customer Service: 1300 753 567 or via Customer Help at target.com.au/help/contact-us

(6.)