


SIZE:44.1X21cm



Massage Vibration Ball

Instruction Manual

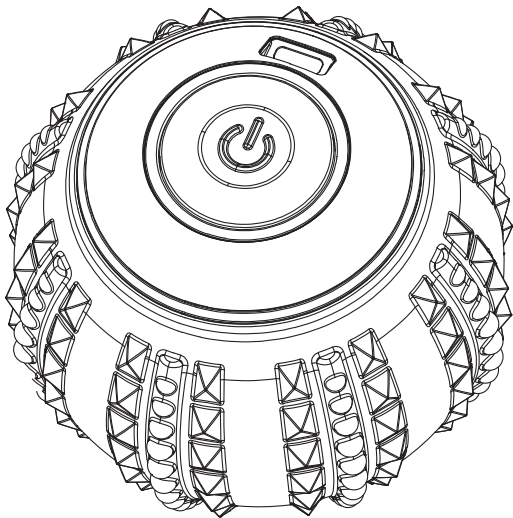

USB Charging Cable

Model H-5916-2

To ensure proper operation of this device, please read this instruction manual carefully before using the device and keep it for further reference.

(1.)

READ AND SAVE IMPORTANT SAFETY INSTRUCTIONS

DANGER- TO REDUCE RISK OF ELECTRIC SHOCK:

1. **DO NOT** reach for an appliance that has fallen into water. Unplug immediately.
2. **DO NOT** use while bathing or in a shower.
3. **DO NOT** place or store appliance where it can fall into water.
4. **DO NOT** place in or drop into water or any other liquid.

WARNING- TO REDUCE THE RISK OF BURNS, FIRE, ELECTRIC SHOCK OR INJURY TO PERSONS:

1. This appliance is not to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction.
2. Children being supervised not to play with the appliance.
3. The appliances must only be charged at safety extra low voltage supply (not provided) corresponding to the marking on the appliance.
4. For the information of charging, please refer the section “Charging” below.
5. This appliance contains batteries that are non-replaceable.
6. Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons. THIS APPLIANCE IS NOT A TOY
7. Use appliance only for its intended use as described in this manual.
8. Never operate this unit if it is not working properly, if it has been dropped or damaged, or dropped into water.
9. **DO NOT** carry this unit by the supply cord or use the cord as a handle.
10. Keep this appliance away from heated surface.
11. **DO NOT** use outdoors.
12. Never drop or insert any object into any opening of this appliance.
13. **DO NOT** operate where aerosol (spray) products are being used or where oxygen being administered.
14. **DO NOT** use over insensitive skin areas or in the presence of poor circulation. The unattended use by children or incapacitated persons may be dangerous.
15. Never use while sleeping and avoid falling asleep while the product is operating.
16. Unplug the adapter after charging and before using the device.
17. **DO NOT** operate continuously for more than 15 minutes. Allow device to rest for 30 minutes before re-using.
18. Keep fingers, hair or any other body part away from the shaft and neck of appliance as pinching may occur.

(2.)

CAUTION- THIS MASSAGER SHOULD NOT BE USED IN THE FOLLOWING CASES:

1. If you have any concerns regarding your health consult your doctor before using this product.
2. Pacemaker recipients should consult a doctor before using this product.
3. Never leave massager unattended, especially if children are present.
4. **DO NOT** fall asleep while using this massager.
5. Massage should be comfortable and pleasant. If pain or discomfort result, discontinue and consult your doctor.
6. Avoid the use of an extension cord because the extension cord may overheat and cause a risk of fire.
7. **DO NOT** use this product before bedtime. The massage could create a stimulating effect and delay sleep.
8. This product should never be used by an individual with a medical condition that would limit the user capacity to operate the controls.
9. Never use this product directly on swollen or inflamed areas of the skin or on areas where skin eruptions are present.
10. **DO NOT** use it in areas close to bones where there is less muscle coverage, such as tibia, dorsum of foot, and back of hand.
11. **DO NOT** use on or around the following sensitive areas: eyes, teeth, breast augmentation, reproductive organs, prostheses, etc.
12. **DO NOT** operate after taking pain killers, sedatives or alcohol. Consult your doctor if under these medications.
13. Consult your doctor before using this product if you are pregnant or have diabetes.
14. **DO NOT** use this product on varicose veins.
15. This appliance is not to be used for medical treatment purposes.

ATTENTION:

1. Fully charged before first use.
2. The massager **cannot be used while charging**.
3. If charging is started during operation, the device will enter **self-protection mode** and automatically shut off.
4. The red indicator light will flash to indicate the device need to be charged. After flash, the light will off, the device is out of power.
5. The red indicator light will flash when charging, once fully charged, the red indicator light will steady on.

CHARGING

- Use the provided **USB** charging cable to connect the massager.
- Use a **DC 5V/2A adapter** (not included) or compatible USB power source.
- **Battery capacity:** 3.7V Lithium-ion battery (1300mAh).
- **Charging time:** approximately **2 hours** for a full charge.
- **Usage time:** around **4 hours** per full charge.
- The massager **cannot operate while charging**.



(3.)

OPERATING INSTRUCTIONS


1. Short press the **Power button**  , the Red indicator light will turn on and display battery power, then go off.
2. The product provides 5 massage intensities + 4 massage modes, you can select by pressing the **Power button**  .

• Press and hold the **Power button**  to turn on the massager;

	Blue light indication		Blue light indication
Gear 1	●○○○○	Mode 1	●●○○○
Gear 2	●●○○○	Mode 2	○●●○○
Gear 3	●●●○○	Mode 3	○○●●○
Gear 4	●●●●○	Mode 4	○○○●●
Gear 5	●●●●●		

- Press again to turn back to **Gear 1**, press the **Power button**  repeatedly to go through different setting.
3. Press and hold the **Power button**  to turn **OFF** when the product is powered on in any mode.
4. The vibration function will automatically shut off after 15 minutes of use. If you need to use it again, please turn it on again. It is recommended to use it for 15 minutes at a time, as not to cause physical discomfort.

Product Description

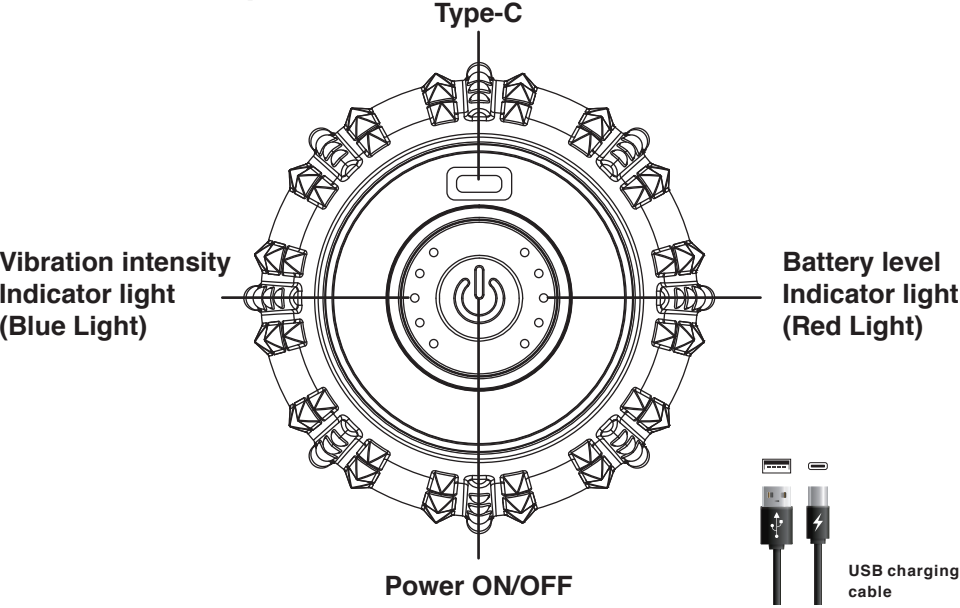

(4.)

CLEANING AND MAINTENANCE

1. Ensure the power switch is turned off and the plug removed before cleaning. Wipe the massager with a **clean, soft cloth**.
2. Do **not** use corrosive cleaning agents such as benzene, thinners, or solvents.
3. Do **not** wipe the control surface with a wet cloth – use a dry cloth only.
4. Do **not** immerse the product in water or any other liquid.
5. If water is accidentally spilled on the product, it may cause electric shock. Keep the product away from **humid areas**. Store it in a **clean, dry place** in its original packaging when not in use. **Keep out of reach of children**.
6. Do **not** dispose of this product with household waste. Follow local **environmental disposal regulations**.

TROUBLE SHOOTING

Problem	Possible Cause / Solution
Massager does not work	Ensure the massager is fully charged.
	Check that the charger and adapter meet the required specifications.
	Confirm the power button has been pressed correctly.
Massager works slowly or weakly	Ensure the massager is fully charged.
	If the contact with your body is loose, adjust the position or pressure for better effect.
Massager suddenly powers off	Check the battery charge level.
	The appliance will automatically power off after 15 minutes of usage.
	If the issue persists, discontinue use and contact customer service.

(5.)

12 Month Warranty

Thank you for your purchase.

Your new product is warrantied to be free from defects in materials and workmanship for the period stated above, from the date of purchase, provided that the product is used in accordance with accompanying recommendations or instructions where provided. This warranty is in addition to your rights under the Australian Consumer Law.

For New Zealand customers, this warranty is in addition to statutory rights observed under New Zealand legislation.

We will provide you with your choice of a refund, repair (where possible) or exchange (availability dependent) for this product if it becomes defective within the warranty period. The business will bear the reasonable expense of claiming the warranty.

This warranty will no longer apply where the defect is a result of alteration, accident, misuse, abuse or neglect.

Please retain your receipt as proof of purchase and contact our Customer Service Centre as listed below, for the entity you purchased this product from, for any difficulties with your product. Warranty claims and claims for expense incurred in returning this product can be addressed to the respective Customer Service Centre.

Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and compensation for any other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure.

Contact for Kmart purchases	Contact for Target purchases
Kmart Australia Ltd. C/- Customer Service Centre L2/ 1 MIDDLE ROAD, CHADSTONE VIC 3148 AUSTRALIA Customer Service: 1800 124 125 (Australia) or 0800 945 995 (New Zealand) or via Customer Help at kmart.com.au	Target Australia Pty Ltd C/- Customer Service Centre L2/ 1 MIDDLE ROAD, CHADSTONE VIC 3148 AUSTRALIA Customer Service: 1300 753 567 or via Customer Help at target.com.au/help/contact-us

(6.)