

anko

Foot Massager Instruction Manual



FMS-125H

Rev 01

To ensure proper operation of this device, please read this instruction manual carefully before using the device and keep it for further reference.

IMPORTANT SAFETY INSTRUCTIONS

THE APPLIANCE IS ONLY TO BE USED WITH THE POWER SUPPLY UNIT PROVIDED WITH THE APPLIANCE.

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

READ ALL INSTRUCTIONS BEFORE USING.

DANGER TO REDUCE THE RISK OF ELECTRIC SHOCK:

- **ALWAYS** unplug this appliance from the electrical outlet immediately after using and before cleaning.

WARNING – TO REDUCE THE RISK OF BURNS, FIRE, ELECTRIC SHOCK, OR INJURY TO PERSONS:

- An appliance should **NEVER** be left unattended
- An appliance should **NEVER** be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- **DO NOT** operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- This product should **NEVER** be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls, or who has sensory deficiencies.
- Use this appliance only for its intended use as described in this manual. **DO NOT** use attachments not provided with the unit.
- **NEVER** operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. contact Consumer Service at the telephone number listed in the Warranty section.
- **DO NOT** carry this appliance by supply cord or use cord as a handle.
- Keep the cord away from heated surfaces.
- **NEVER** operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- **NEVER** drop or insert any object into any opening.
- **DO NOT** use outdoors.
- **DO NOT** operate where aerosol (spray) products are being used or

IMPORTANT SAFETY INSTRUCTIONS

where oxygen is being administered.

- To disconnect, turn all controls to the off position, then remove plug from outlet.
- Use heated surfaces carefully. May cause serious burns. **DO NOT** use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- **DO NOT** use massager in close proximity to loose clothing or jewelry.
- Keep long hair away from massager while in use.
- **DO NOT** stand on or in appliance. Use only while seated.
- **CAUTION:** All servicing of this massager must be performed by authorised personnel only. For service, contact Consumer Service at the telephone number listed in the Warranty section.

SAVE THESE INSTRUCTIONS

CAUTION

PLEASE READ ALL INSTRUCTIONS CAREFULLY BEFORE OPERATING.

- **WARNING** Only use the Power Supply that was provided with this appliance. Model No. XH2400-1000ENCs.
- This product is not intended for medical use. It is intended only to provide a luxurious massage.
- Consult your doctor prior to using this product, if:
 - You are pregnant
 - You have a pacemaker
 - You have any concerns regarding your health
- **NOT** recommended for use by diabetics.
- **NEVER** leave the appliance unattended, especially if children are present.
- **NEVER** cover the unit when it is in operation.
- **DO NOT** use this product for more than 15 minutes at a time.
- Extensive use could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.
- **NEVER** use this product directly on swollen or inflamed areas or skin eruptions or where sores are present.
- **DO NOT** use this product as a substitute for medical attention.
- **DO NOT** use this product before bed. The massager has a stimulating effect and can delay sleep.
- **NEVER** use this product while in bed.
- This product should **NEVER** be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls.
- This unit should **NOT** be used by children or invalids without adult supervision.
- This appliance is intended for household use only.
- This product is not for use in automobiles.
- **NOTE:** This unit may turn off after an extended period of use. If your unit turns off (LEDs may remain on), unplug the unit and let it rest for 30 minutes before plugging in again and continuing use. Also, the heating function and LEDs may turn off after a period of use due to the


YOUR FOOT MASSAGER

DEEP-KNEADING SHIATSU MESSAGE NODES

Provide an invigorating massage

SOOTHING HEAT

For a relaxing massage

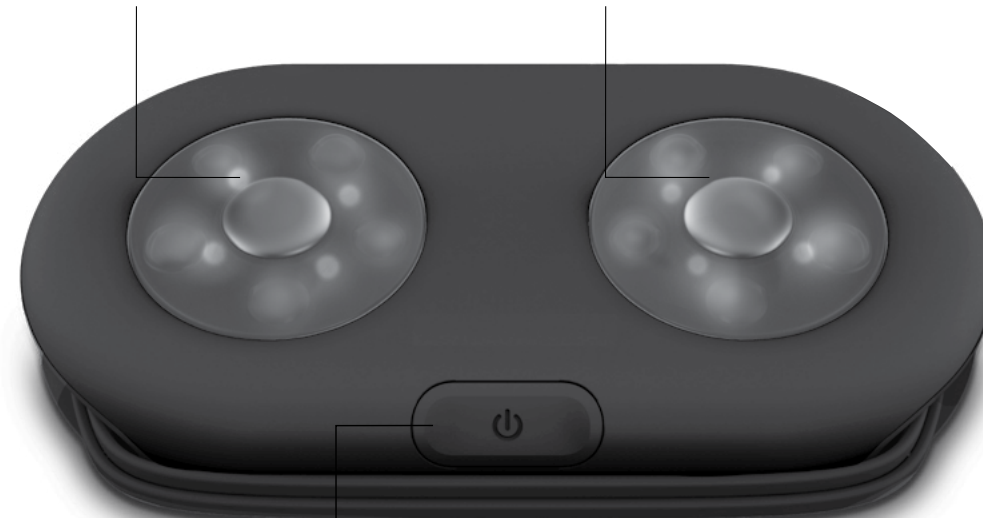

TOE-TOUCH CONTROL

Activates massage and heat

POWER ADAPTOR

Power adaptor included



INSTRUCTIONS FOR USE

Fig. 1

1. Connect adaptor to massager and then plug into a power outlet.
2. It is highly recommended that the foot massager be used with socked feet.
DO NOT use the unit with shoes.



Fig. 1

3. Place massager on the floor in front of your chair or seat.
4. To activate massage action, press the power button (Fig. 1).
Tip: By changing the pressure of your feet on the massager, it will change the feel of the massage. If the unit tickles you, then you may need to apply more pressure or you may need to become accustomed to the sensation. If the massage feels too strong, you may need to lighten the pressure of your feet on the unit.
5. To activate soothing heat with the massage, press the power button a second time. Allow a few moments for the heating action to warm up.
6. When finished, press a third time and the unit will turn off.
Unplug the unit.

NOTE: As a safety feature, the soothing heat function cannot be used alone.

CAUTION: To avoid injury and irreparable damage, always remain seated when using massager. **NEVER** stand on the unit.

MAINTENANCE

TO STORE

Place massager in its box or in a cool, dry place. Avoid contact with sharp edges or pointed objects. To avoid breakage, **DO NOT** wrap the power cord over the unit, use the cord storage around the massager.

DO NOT hang the unit by the cord.

TO CLEAN

Unplug the unit and allow it to cool before cleaning. Clean only with a soft, slightly damp cloth. **NEVER** allow water or any other liquids to come into contact with the unit.

- **DO NOT** immerse in any liquid to clean.
- **NEVER** use abrasive cleaners, brushes, or strong chemicals that may or may not be flammable and/or damaging to the product to clean.
- **DO NOT** attempt to repair the massager. There are no user-serviceable parts.
- For service, contact Consumer Service at the telephone number listed in the Warranty section.

12 Month Warranty

Thank you for your purchase.

Your new product is warrantied to be free from defects in materials and workmanship for the period stated above, from the date of purchase, provided that the product is used in accordance with accompanying recommendations or instructions where provided. This warranty is in addition to your rights under the Australian Consumer Law.

For New Zealand customers, this warranty is in addition to statutory rights observed under New Zealand legislation.

We will provide you with your choice of a refund, repair (where possible) or exchange (availability dependent) for this product if it becomes defective within the warranty period. The business will bear the reasonable expense of claiming the warranty.

This warranty will no longer apply where the defect is a result of alteration, accident, misuse, abuse or neglect.

Please retain your receipt as proof of purchase and contact our Customer Service Centre as listed below, for the entity you purchased this product from, for any difficulties with your product. Warranty claims and claims for expense incurred in returning this product can be addressed to the respective Customer Service Centre.

Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and compensation for any other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure.

Contact for Kmart purchases	Contact for Target purchases
Kmart Australia Ltd. C/- Customer Service Centre 690 Springvale Road, Mulgrave Vic 3170 Customer Service: 1800 124 125 (Australia) or 0800 945 995 (New Zealand) or via Customer Help at kmart.com.au	Target Australia Pty Ltd C/- Customer Service Centre 2 Kendall Street, Williams Landing, Vic, 3027 Customer Service: 1300 753 567 or via Customer Help at target.com.au/help/contact-us