



Compression Calf Instruction Manual



YY-A09C

Rev 01

To ensure proper operation of this device, please read this instruction manual carefully before using the device and keep it for further reference.

IMPORTANT SAFETY INSTRUCTIONS - RETAIN FOR FUTURE USE

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING.

DANGER: TO REDUCE THE RISK OF ELECTRIC SHOCK:

- Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
- **DO NOT** reach for a unit that has fallen into water. Unplug it immediately.

WARNING: TO REDUCE THE RISK OF BURNS, FIRE, ELECTRIC SHOCK, OR INJURY TO PERSONS:

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- An appliance should **NEVER** be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts or attachments.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- **DO NOT** use outdoors.
- Use this appliance only for its intended use as described in this manual. **DO NOT** use attachments any attachments not provided with the unit.
- **DO NOT** carry this appliance by the supply cord or use cord as a handle.
- **NEVER** operate this appliance if it has a damaged USB cable, if it is not working properly, if it has been dropped or damaged, or if it has been dropped into water. Return the appliance for examination and repair.
- Keep the cord away from heated surfaces.
- **NEVER** operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- **NEVER** drop or insert any object into any opening.
- **DO NOT** operate where aerosol (spray) products are being used or where oxygen is being administered.
- **DO NOT** operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- **DO NOT** stand on or in appliance. Use only while seated.
- Keep long hair away from massager while in use.

SAVE THESE INSTRUCTIONS

- **DO NOT** use massager in close proximity to loose clothing or jewellery.
- **NEVER** operate on a soft surface such as a bed or couch where the air openings may be blocked.
- Prevent unintentional starting. Ensure the switch is in the off position on the control box before picking up or carrying the appliance. **DO NOT** carry the unit with your fingers on the control box buttons.
- Turn off the control box from the appliance before making any adjustments, changing accessories, or storing appliance. Such preventive safety measures reduce the risk of starting the appliance accidentally.
- Recharge only with the USB cable specified by the manufacturer. A USB cable that is suitable for one type of electrical device may create a risk of fire when used with another.
- Use appliances only with specifically designated control box. Use of any other may create a risk of injury and fire.
- When massager is not in use, keep it away from other metal objects, like paper clips, coins, keys, nails, screws or other small metal objects, that can make a connection from one terminal to another. Shorting the battery terminals together may cause burns or a fire.
- Under abusive conditions, liquid may be ejected from the battery; avoid contact. If contact accidentally occurs, flush with water. If liquid contacts eyes, seek medical help.
- Liquid ejected from the battery may cause irritation or burns.
- **DO NOT** use the control box or appliance that is damaged or modified. Damaged or modified devices may exhibit unpredictable behavior resulting in fire, explosion or risk of injury.
- **DO NOT** expose the control box or appliance to fire or excessive temperature.
Exposure to fire or temperature over 130°C may cause explosion.
- Follow all charging instructions and **DO NOT** charge the battery or appliance outside of the temperature range specified in the instructions. Charging improperly or at temperatures outside of the specified range may damage the battery and increase the risk of fire. Product operating and charging range: 0°C – 40°C
- This appliance contains batteries that are only replaceable by skilled persons.
- This appliance contains batteries that are non-replaceable by the user.



SAVE THESE INSTRUCTIONS

CAUTION: PLEASE READ ALL INSTRUCTIONS CAREFULLY BEFORE OPERATING.

- Consult your doctor prior to using this product, if:
 - You are pregnant
 - You have a pacemaker
 - You have any concerns regarding your health
- **NOT** recommended for use by diabetics.
- **NEVER** leave the unit unattended, especially if children are present.
- **NEVER** cover the unit when it is in operation.
- **DO NOT** use this product for more than 15 minutes at a time.
- Extensive use could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow unit to cool before operating.
- **NEVER** use this product directly on swollen or inflamed areas or skin eruptions.
- **DO NOT** use this product as a substitute for medical attention.
- **DO NOT** use this product before bed. The massager has a stimulating effect and can delay sleep.
- This product should **NEVER** be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls.
- **NOT** for use by children.
- This appliance is intended for household use only.
- This product should **NEVER** be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls or who has sensory deficiencies.
- **NOTE:** This unit will shut off automatically after 15 minutes of use. If your unit overheats, turn off the unit and let cool for 30 minutes before turning on again and continuing use.

INSTRUCTIONS FOR USE

- 1. POWER SWITCH**
Long press to power on/off.
- 2. BATTERY INDICATOR**
Green: 100%
Blue: 60%
Red: 20%
- 3. MODE INDICATOR**
Four Massage Modes:
A: Alternating
b: Full
c: Sequential
d: Pulsing
- 4. PRESSURE INDICATOR**
5 levels of pressure:
20% / 40% / 60% / 80% / 100%
- 5. DECREASE PRESSURE**
Short press the '-' button to decrease the compression pressure level.
- 6. INCREASE PRESSURE**
Short press the '+' button to increase the compression pressure level.
- 7. SHORT PRESS TO SWITCH MODES.**
It has A/B/C/D, 4 modes to adjust.
Automatic operation of default Mode A after 2 seconds.
- 8. CHARGING PORT**
Type-C charging port, input 5V-1A

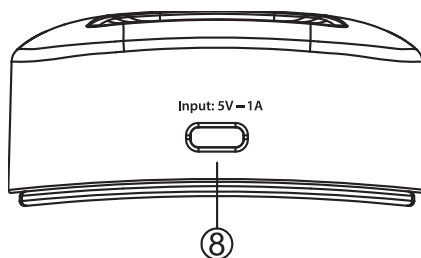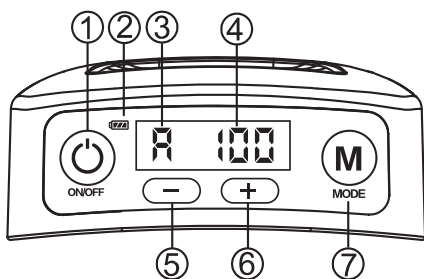


CHARGING THE COMPRESSION UNIT

When the battery is low, the battery symbol will flash red.

1. To charge the product, plug the charging cable into a USB socket and connect the cable to the USB-C power outlet source.
2. Once the charging cable is connected, the battery indicator will flash green when the battery is charging.


NOTE It might take approximately 3.5 hours to fully charge the battery. The battery run time will depend on what pressure level and massage mode has been used.

INSTRUCTIONS FOR USE

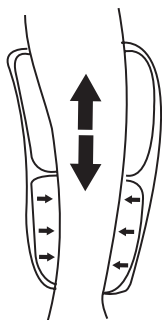


MODE A: ALTERNATING

The first airbag is inflated to the desired pressure and maintains the pressure for 5 seconds. The air is then released for 2 seconds. The second airbag begins to inflate to the desired pressure and maintains the desired pressure for 5 seconds. The air is then released for 2 seconds. Both airbags inflate simultaneously and maintain pressure for 5 seconds. The pattern continues to repeat until auto shutoff at 15 minutes.

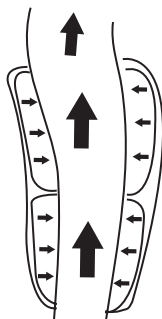


MODE B: FULL

The two airbags are inflated to the desired pressure at the same time. The airbags hold the pressure for 5 seconds and then release the air for 3 seconds. This pattern continues to repeat until auto shutoff at 15 minutes.

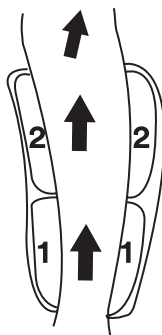

MODE C: SEQUENTIAL

The lower airbag is inflated to the desired pressure, and the pressure is held for 3 seconds. The upper airbag then inflates to the desired pressure. Air pressure is held in both airbags for 5 seconds and then released for 5 seconds. This pattern continues to repeat until auto shutoff at 15 minutes.

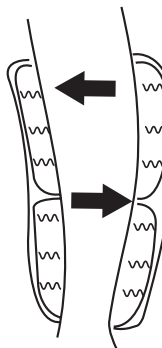

MODE D: PULSING

The two airbags are inflated to the desired pressure at the same time, then the airbags immediately release the air for 1 second. The airbags then continue to inflate to the desired pressure. This pattern continues to repeat until auto shutoff at 15 minutes.

MAINTENANCE

TO CLEAN

- Switch off the massager before cleaning.
- Clean only with a soft, slightly damp sponge.
- **NEVER** allow water or any other liquids to come into contact with the appliance.
- **DO NOT** immerse in any liquid to clean.
- **NEVER** use abrasive cleaners, brushes, glass / furniture polish, paint thinner etc to clean.

TO STORE

- Ensure the product is turned off and then place in its box or in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface. To avoid breakage.
- **DO NOT** wrap the USB cable around the appliance.
- **DO NOT** hang the unit by the USB cable.

BATTERY DISPOSAL



At the end of the appliance's life, it must be disposed of correctly.

- Contact your local waste management authority for more information on Battery Recycling and disposal regulations in your area.
- Do not dispose of electrical appliances or batteries as unsorted municipal waste; use separate collection facilities.
- Contact your local council for information regarding the collection systems available. If electrical appliances are disposed of in landfills or dumps, hazardous substances can leak into the groundwater and get into the food chain, damaging your health and well-being.
- Please contact your local council or regional authority for more information on the collection, reuse, and recycling programs.

12 Month Warranty

Thank you for your purchase.

Your new product is warranted to be free from defects in materials and workmanship for the period stated above, from the date of purchase, provided that the product is used in accordance with accompanying recommendations or instructions where provided. This warranty is in addition to your rights under the Australian Consumer Law.

For New Zealand customers, this warranty is in addition to statutory rights observed under New Zealand legislation.

We will provide you with your choice of a refund, repair (where possible) or exchange (availability dependent) for this product if it becomes defective within the warranty period. The business will bear the reasonable expense of claiming the warranty.

This warranty will no longer apply where the defect is a result of alteration, accident, misuse, abuse or neglect.

Please retain your receipt as proof of purchase and contact our Customer Service Centre as listed below, for the entity you purchased this product from, for any difficulties with your product. Warranty claims and claims for expense incurred in returning this product can be addressed to the respective Customer Service Centre.

Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and compensation for any other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure.

Contact for Kmart purchases	Contact for Target purchases
Kmart Australia Ltd. C/- Customer Service Centre 690 Springvale Road, Mulgrave Vic 3170 Customer Service: 1800 124 125 (Australia) or 0800 945 995 (New Zealand) or via Customer Help at kmart.com.au	Target Australia Pty Ltd C/- Customer Service Centre 2 Kendall Street, Williams Landing, Vic, 3027 Customer Service: 1300 753 567 or via Customer Help at target.com.au/help/contact-us