




Nasal Aspirator Instruction Manual

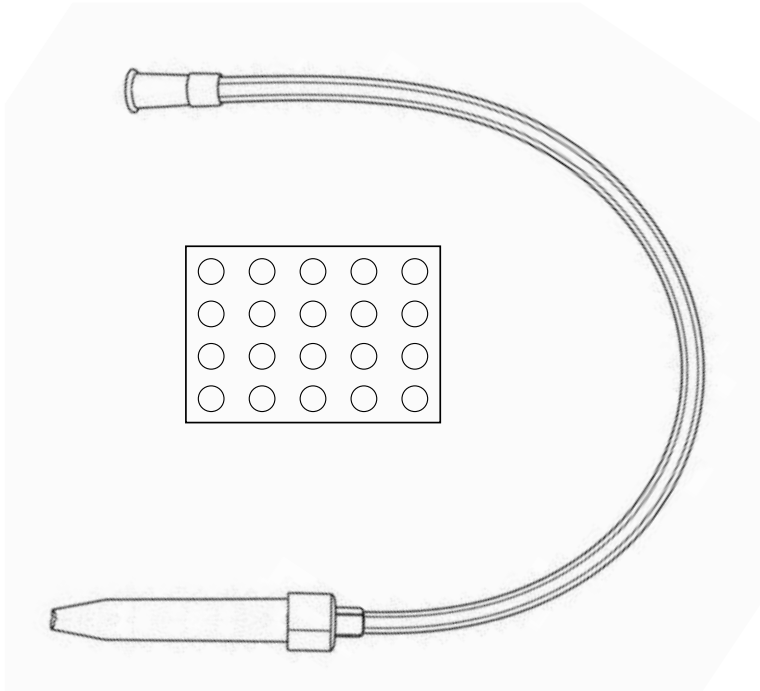

IMPORTANT! READ CAREFULLY AND KEEP FOR FUTURE REFERENCE

Read all instructions **BEFORE** assembly and **USE** of product

Note: Drawings in the manual are for illustration purposes only. Actual product may vary slightly.
The manufacturer reserves the right to change any specification or features without prior notice.

1. SCOPE OF APPLICATION

THIS UNIT IS A MANUALLY OPERATED NASAL ASPIRATOR DESIGNED FOR DOMESTIC USE FOR NEW BORN AND ABOVE.
TO SAFELY AND GENTLY HELP REMOVE SECRETIONS AND MUCUS.

2. PARTS DESCRIPTION

PARTS	FUNCTION	QTY
LARGE TUBE	FITS THE BABY'S NOSTRILS AND SAFELY GUIDES NASAL SECRETIONS	1
LONG TUBE	CONNECTS TO THE LARGE TUBE AND PREVENTS BACKFLOW OF NASAL SECRETIONS	1
FILTER SPONGE	FILTERS AND BLOCKS NASAL SECRETIONS TO PREVENT LIQUID FROM FLOWING BACK TO THE SUCTION NOZZLE	20
INSTRUCTION MANUAL	GUIDES PRODUCT ASSEMBLY, USE, CLEANING, DISINFECTION AND SAFETY PRECAUTIONS TO ENSURE STANDARD OPERATION	1

3. HOW TO USE

PREPARATION

- › BEFORE FIRST USE AND AFTER EACH USE, DISASSEMBLE ALL PARTS, DISINFECT THEM, AND AIR-DRY. INSERT THE SPONGE INTO THE LONG TUBE, THEN CONNECT THE LARGE TUBE TO THE LONG TUBE. MAKE SURE ALL CONNECTIONS ARE TIGHTLY FITTED.
- › SECRETIONS: LAY THE BABY FLAT WITH THE HEAD SLIGHTLY ELEVATED. PUT 1-2 DROPS OF INFANT-SPECIFIC SALINE INTO EACH NOSTRIL. GENTLY MASSAGE THE NASAL ALAE, AND WAIT 30 SECONDS TO 1 MINUTE TO DILUTE STICKY SNOT AND SOFTEN DRY BOOGERS.FOR DRY AND HARD BOOGERS, APPLY A WARM TOWEL TO THE NOSE FOR 1-2 MINUTES TO SOFTEN THEM.
- › PLACE THE BABY IN SUPINE / SEMI-SUPINE POSITION WITH THE HEAD ELEVATED ABOUT 30°. THE PARENT FIXES THE BABY'S HEAD TO AVOID SHAKING.

OPERATION STEPS

- › THE PARENT GENTLY SUCKS THE SUCTION NOZZLE FIRST TO FEEL THE SUCTION STRENGTH AND ENSURE IT IS MILD AND CONTROLLABLE.
- › GENTLY PLACE THE LARGE TUBE AT THE ENTRANCE OF THE NOSTRIL, ALIGN WITH THE NASAL PASSAGE, AND KEEP THE LARGE TUBE CLOSE TO THE NOSTRIL. DO NOT INSERT DEEP TO TOUCH THE NASAL MUCOSA.
- › SUCK SLOWLY AND GENTLY WITH THE SUCTION NOZZLE UNTIL A SLIGHT SUCTION SOUND IS HEARD. SINGLE SUCTION ≤ 3 SECONDS. SECRETIONS WILL FLOW INTO THE LARGE TUBE.
- › AFTER ONE SIDE IS CLEANED, PULL OUT THE LARGE TUBE AND REPEAT THE OPERATION ON THE OTHER NOSTRIL.
- › CLEAN RESIDUE GENTLY WIPE AROUND THE NOSTRILS WITH SOFT TISSUE AFTER SUCTION. DO NOT RUB HARD.

TO CLEAN

- › DISASSEMBLE ALL PARTS. POUR OUT SECRETIONS FROM THE LARGE TUBE, AND DISCARD THE SPONGE.CLEAN ALL PARTS WITH WARM WATER AND INFANT-SPECIFIC CLEANING AGENT. GENTLY BRUSH THE INNER WALL OF THE LARGE TUBE WITH A SOFT BRUSH.RINSE WITH CLEAN WATER, AIR-DRY IN A VENTILATED PLACE, ASSEMBLE AND STORE IN A DRY AND CLEAN STORAGE BOX TO AVOID DUST CONTAMINATION.

WARNING

1. EXCESSIVE SUCTION, IMPROPER ASSEMBLY, DAMAGED OR AGED PARTS, OR INADEQUATE CLEANING MAY CAUSE BACKFLOW OF SECRETIONS, POSING A RISK OF INFECTION TO THE OPERATOR.
2. FOR INFANT NASAL USE ONLY. DO NOT USE ON EARS, MOUTH OR OTHER BODY PARTS.
3. DO NOT FORCE THE LARGE TUBE INTO NOSTRILS OR INSERT DEEP INTO THE NASAL CAVITY TO AVOID INJURY AND BLEEDING OF NASAL MUCOSA.
4. SUCTION MUST BE GENTLE. DO NOT SUCK HARD OR CONTINUOUSLY FOR A LONG TIME. USE ≤ 3 TIMES A DAY. EXCESSIVE STIMULATION MAY CAUSE EDEMA OF NASAL MUCOSA.
5. STOP IMMEDIATELY IF THE BABY CRIES OR RESISTS. DO NOT OPERATE FORCIBLY.
6. FOR PERSONAL USE ONLY. DO NOT SHARE WITH OTHERS TO AVOID CROSS-INFECTION.
7. CLEAN AND DISINFECT AFTER EACH USE TO PREVENT BACTERIAL GROWTH.
8. STORE IN A DRY AND COOL PLACE OUT OF REACH OF INFANTS. DO NOT EXPOSE TO SUNLIGHT OR HIGH TEMPERATURE BAKING.
9. DO NOT USE IF THE NASAL CAVITY HAS WOUNDS, BLEEDING, ULCERS OR ACUTE INFLAMMATION.
10. STOP USING AND CONSULT A DOCTOR IMMEDIATELY IF SECRETIONS ARE BLOODY, YELLOW-GREEN PURULENT, OR ACCOMPANIED BY FEVER, DIFFICULTY BREATHING, ETC.
11. DO NOT ALLOW INFANTS TO USE BY THEMSELVES TO AVOID SWALLOWING PARTS OR DAMAGING THE NASAL CAVITY.

KEYCODE:

K: 43-728-435 | T: 72-511-183

MADE IN CHINA

CUSTOMER SERVICE:

KMART AU: 1800 124 125

KMART NZ: 0800 945 995

TARGET AU: 1300 753 567