


FRONT COVER
7cm



anko

Lactation Massager Instruction Manual

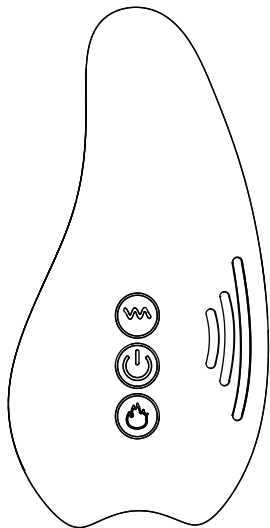

THE ILLUSTRATIONS SHOWN IN THIS MANUAL ARE FOR REFERENCE ONLY. THE MANUFACTURER RESERVES THE RIGHT TO CHANGE ANY SPECIFICATION OR FEATURE WITHOUT FURTHER NOTICE. IMPORTANT! PLEASE READ CAREFULLY BEFORE USE AND KEEP FOR FUTURE REFERENCE.

Intended Use:

This product is designed for external massage of the breast area to provide comfort during breastfeeding.

- › For adult use only
- › For external use only
- › This product is not a medical device and is not intended to treat, cure, or prevent any medical condition. If you experience pain, mastitis, or breastfeeding concerns, consult a qualified healthcare professional.

Warnings:

General Safety

- › Do not use while sleeping.
- › Do not apply excessive pressure during use.
- › This product is not intended for use by children.
- › Children should be supervised to ensure they do not play with the appliance.
- › Do not use if the device is damaged or malfunctioning.

Heat Safety

- › The device has a heated surface. Use caution if sensitive to heat.
- › Do not use on broken, irritated, or sensitive skin.
- › Do not use heat continuously for extended periods.
- › Check skin condition regularly during use.
- › Discontinue use immediately if discomfort, redness, or burning sensation occurs.

Battery & Electrical Safety

- › This product contains a non-replaceable rechargeable lithium-ion battery.

- › Do not disassemble, crush, puncture, or modify the device or battery.
- › Do not expose to fire, high temperatures, or direct sunlight (e.g. inside a car).
- › Do not use or charge if the device is swollen, leaking, overheating, or emitting unusual odours.
- › Do not charge the device unattended.
- › Ensure the device is switched off before charging.

Water & Environment

- › Do not immerse the device in water.
- › Do not use or charge near water or in damp environments.
- › Keep away from liquids, fire, and extreme heat. Consult a healthcare professional before use if you have any medical conditions.


Package contents:

- 1 × Lactation Massager
- 1 × Instruction Manual
- 1 × USB-A to DC Charging Cable

Product specifications:

- Input: 5V= 5W
- Battery Capacity: 600mAh, 3.7V rechargeable
- Lithium-ion battery (non-replaceable)
- Water resistance rating: IPX6
- Vibration Modes: 10 selectable vibration patterns across 3 modes (Constant, Pulse, and Oscillating).
- Heat settings: 3 levels – Low, Medium, High

Product diagram:

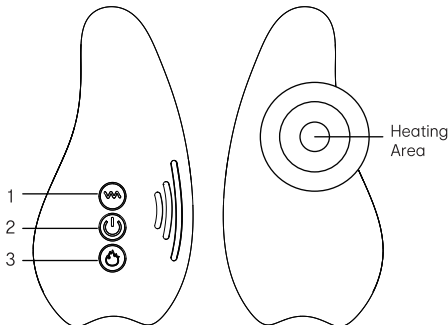

1. Vibration Mode Button
2. Power Button
3. Heat Button

Instructions for use:

Power On/Off

- › Long press the power button for 2 seconds to turn the device on or off.
- › Turning the device off will stop all functions, including vibration and heat.

Vibration Function

- › Short press the Vibration Mode Button to cycle through the available vibration modes and intensity levels.
- › To stop vibration, long press the Power Button for 2 seconds to turn the device off.

Heat function

The heat levels cycle in sequence.

Short press the Heat Button to cycle through the heat settings:

- › First press: Low heat
- › Second press: Medium heat (LED flashes twice)
- › Third press: High heat (LED flashes three times)
- › Fourth press: Heat off.

How to charge:

- › Use the supplied USB-A to DC charging cable with a suitable 5V DC 1A power adaptor (not included) to charge the device.
- › The indicator light will flash while charging.
- › The indicator light will remain on when fully charged.
- › Fully charge before first use.

Cleaning & Maintenance:

- › Clean before and after each use.
- › Wipe the surface with a soft damp cloth.
- › Use mild detergent only if required.
- › Do not use harsh chemicals or abrasive cleaners.
- › Do not immerse in water. Ensure the device is completely dry before storage.
- › For personal use only – do not share.
- › Ensure the device is clean before skin contact.
- › Do not use on infected, damaged, or irritated skin.
- › Store in a cool, dry place.
- › Avoid direct sunlight and high temperatures.
- › Do not place heavy objects on the device.

Notes: This product is intended for personal comfort use only and is not a medical device. If you experience persistent pain, mastitis, or breastfeeding difficulties, consult a healthcare professional.

12 Months Warranty

Thank you for your purchase.

Your new product is warrantied to be free from defects in materials and workmanship for the period stated above, from the date of purchase, provided that the product is used in accordance with accompanying recommendations or instructions where provided. This warranty is in addition to your rights under the Australian Consumer Law.

For New Zealand customers, this warranty is in addition to statutory rights observed under New Zealand legislation.

We will provide you with your choice of a refund, repair (where possible) or exchange (availability dependent) for this product if it becomes defective within the warranty period. The business will bear the reasonable expense of claiming the warranty.

This warranty will no longer apply where the defect is a result of alteration, accident, misuse, abuse or neglect.

Please retain your receipt as proof of purchase and contact our Customer Service Centre on 1800 124 125 (Australia) or 0800 945 995 (New Zealand) or alternatively, via the Customer Service number below for the entity from which you purchased the product for any difficulties with your product. Warranty claims and claims for expenses incurred in returning this product can be addressed to our Customer Service Centre at L2/1 MIDDLE ROAD, CHADSTONE VIC 3148.

Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and compensation for any other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure.

Model: JLR-82114

KEYCODE: K: 43-728-688 | T: 72-511-589

MADE IN CHINA

CUSTOMER SERVICE:
KMART AU: 1800 124 125
KMART NZ: 0800 945 995
TARGET AU: 1300 753 567

